



Low-Dead-Space Injection System

Proma inject Syringe + ZELOSTAT™ 35G needle

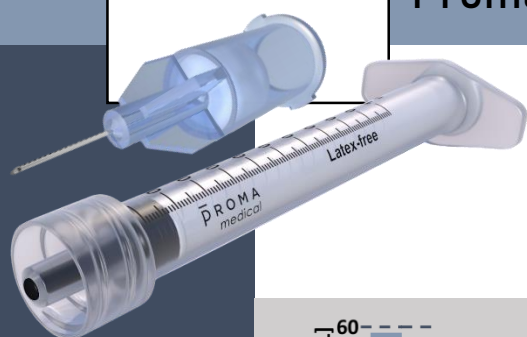


Low-Dead-Space (LDS) Luer-Lock 1mL Proma inject Syringe AND 35G 4mm Low-Dead-Space Hypodermic Needle by ZELOSTAT™ for Botulinum Toxin Micro Injections

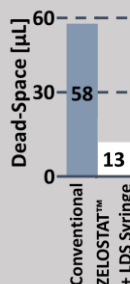

More Product for Your Patient

The Low-Dead-Space technology based on a soft polymer ring filling any extra space, dramatically reducing product waste. From the first use you can already save up to 45µL botulinum toxin solution. Changing the needles 3 times will already save over 100µL of product or one extra injection point you can provide your patient.

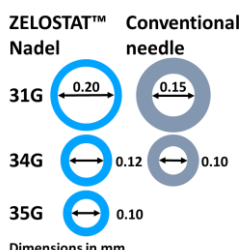

Keep Your Routine

With an inner diameter of 0.1mm, our needle offers a flow and plunger pressure comparable to a conventional 34G, allowing you to maintain your current injection technique. The injector feels the same, but the patient experiences less needle.



Safety and Comfort

The luer-lock adapter tightly secures the injection needle on the syringe allowing the injector to work with a wider range of pressure on the plunger without risking the needle to detach. The specially designed plunger slides without blocking and ensures high level of control during the injection. And, our 35G needle is approximately 25% thinner than conventional 34G needles and can therefore reduce the risk of injection-related side effects, thereby increasing patient comfort during the injection.

Innovation That Pays

The efficiency in our state-of-the-art manufacturing facilities allows us to cut costs dramatically while maintaining the highest standards set by the American FDA and European EMA. We are proud to be able to forward the cost efficiency of our innovations to our customers offering the Proma inject Syringes at a price significantly below market average.

Visit us and order here:



The information provided in this data sheet is for general guidance only and may be subject to change. It is the user's responsibility to ensure the product is suitable for its intended purpose and used according to the manufacturer's instructions. Single use only. Do not reuse or re-sterilize.

PromaInject is part of PromaMedical GmbH
Wernher-von-Braun-Straße 10
85640 Putzbrunn – Munich,
Germany

T: +49 (0)89 7474713-0
F: +49 (0)89 7474713-1
M: info@promainject.com
W: www.promainject.com

e-mail us:

